ALL-IN-ONE SONIC WATER JET

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CONTINENTAL CONAIR LIMITED CHINA CREATIVE DEPARTMENT

Date	6/5/18
Customer	CONAIR-TW
Model No.	CSWJ1W
Description	IB
Job No.	16CFO05
Reference No.	IB-16/270A
Revision No.	1
Color	1C+10(Black)
Open Size	920mm(W) X 214mm(H)
Folded Size	92mm(W) X 214mm(H)
Material	80gsm Artpaper
Prepared by	Penny Peng
Approved by	Jimmy Ng
Vendor	Gimell
Supplier	

折疊方法
Folding method:
20 pages (10 x 2)

Warranty Card

Dear Consumer,
Thank you for your purchase of this products. This warranty card is an important proof for warranty service. Please complete it earnestly. Thanks!

We warranty that under normal use, if the quality or parts of your product should prove to be defective within 1 year from your purchase date, our Consumer Service Center will repair or replace it at our option, without charge to you. To obtain warranty service, please complete and send the warranty card and purchase proof copies to the Service Center listed in the Instruction Book **WITHIN ONE MONTH**.

Please fill in the following information:

Name: Mr./Miss/Ms. _____ Zip Code : _____

Address: _____

Tel: _____ Email: _____

Product Name: _____

Model No.: _____


Store of Purchase: _____

Purchase Date: _____

Warranty Terms:

- Your product has been manufactured to the strictest specifications and has been designed for use with authorized accessories and replacement parts. This warranty expressly excludes any defects or damages caused by accessories, replacement parts or repair service other than those which have been authorized by us.
- This warranty does not cover any damage caused by accident, misuse, shipment or other than ordinary household use.
- This warranty excludes all incidental and consequential damages.
- In addition to return of warranty card, the consumer has to maintain the original proof of purchase in order to obtain the warranty benefits. In the event that you do not have proof of purchase date, the purchase date for purposes of this warranty will be the date of manufacture.


INTERPLAK® ALL-IN-ONE SONIC WATER JET SYSTEM DIAGRAM

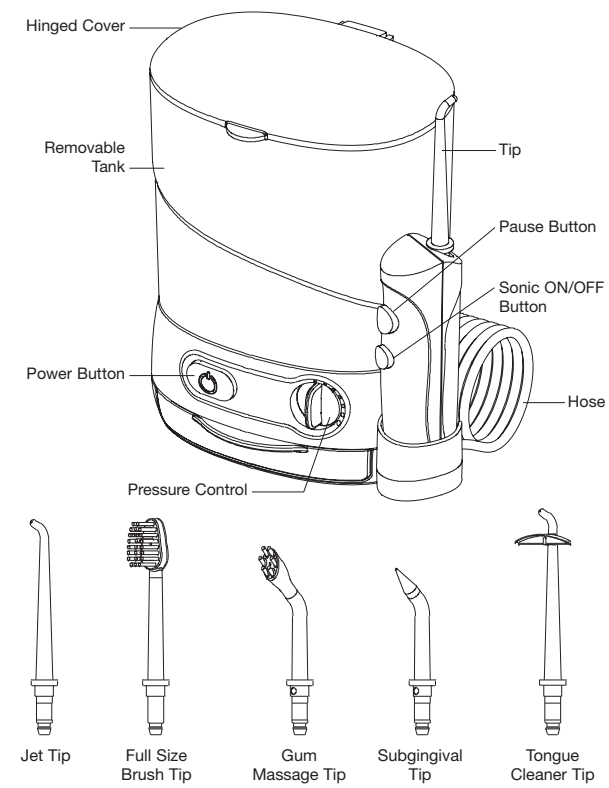GETTING STARTED BATTERY INSTALLATION

Follow these instructions to install battery.

Use a coin to turn bottom cap counterclockwise to Open position; remove cap.

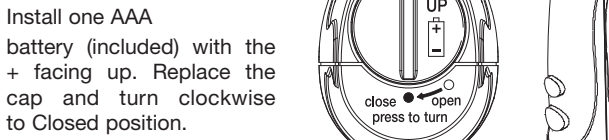FILLING THE WATER TANK

Remove the water tank from the base and fill with lukewarm water. Put the tank back on base and press down firmly to seat it.

INSERTING AND REMOVING TIPS

Select the desired tip and push it firmly into place on the handle. To remove, simply pull the tip off.

PRIMING THE WATER JET

On first use, or if the unit has not been used for a while, the water jet needs to be primed with water.

Fill the water tank to the **max** fill line (900ml) and turn the pressure dial to the maximum setting.

Point the water jet attachment towards the sink and turn the unit on. Within 1 minute, water will begin to flow through the attachment.


Reduce the pressure to the lowest setting (or your desired level) before use.

USING YOUR INTERPLAK® ALL-IN-ONE SONIC WATER JET SYSTEM

- Adjust the Pressure Control knob to the lowest setting (small dot) for first-time use. With continued use, gradually increase pressure to a comfortable setting, or as instructed by your dentist.
- Press the Sonic ON/OFF button on the handle (skip this step if you don't want to use the sonic feature).
- Aim the tip into the sink and press the Power button on the base to start the unit.
- When water starts to flow, press the Pause button and place tip in mouth following the technique described on page 8.
- Lean over the sink and release the Pause button. Press Pause again to stop the flow of water.
- When finished, press the Power button to shut off the unit, and press the Sonic ON/OFF button to turn off the sonic function.

RECOMMENDED TECHNIQUE – WATER JET

Slightly close lips as shown to avoid splashing.

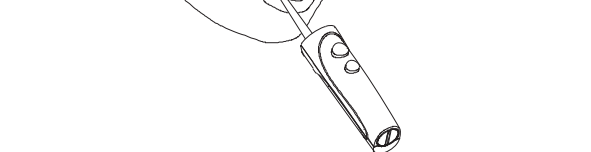Point the water stream at a 90° angle to your gum line.

Slowly move the tip along the gum line.

Follow the recommended pattern until all the teeth have been cleaned.

MASSAGING TIP

Use the massage tip with the sonic vibration function to gently massage and stimulate gums. Press the massaging nodes gently against the gum and follow the method illustrated above. The tip can be used with or without water stream.

When finished, turn off the unit by pressing the Power button on the base and empty the water tank. Be sure to also turn off the Sonic function by pressing the Sonic ON/OFF button on the handle.

RECOMMENDED TECHNIQUE – SONIC BRUSHING



Apply toothpaste. Insert brush tip into mouth, press Sonic button and brush your teeth as you normally do.

RINSING USING WATER STREAM




Press and hold the Pause button and place brush in mouth, slightly closing lips to avoid splashing. Press the Power button, release the Pause button, and rinse along the gum line as shown. Slightly close lips as shown to avoid splashing.

Release the Pause button and rinse as shown along the gum line.









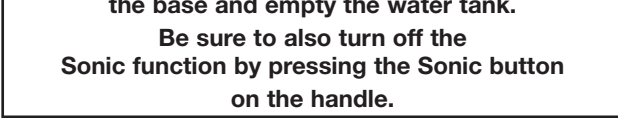
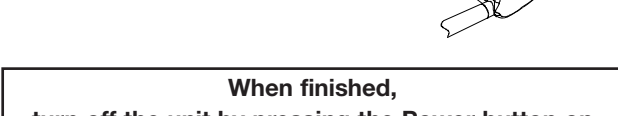
ORTHODONTIC CLEANING

Using the water jet tip, follow the recommended technique described on page 7 to clean around orthodontic brackets one tooth at a time.

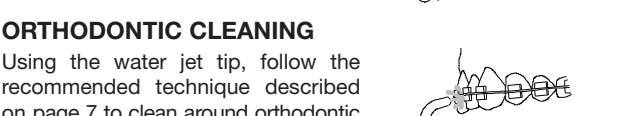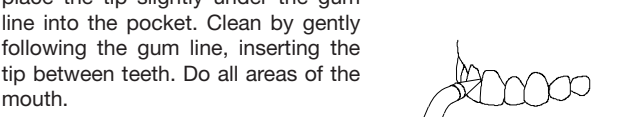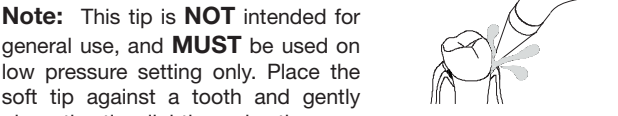USING MOUTHWASH AND OTHER SOLUTIONS





Your Interplak System can be used with mouthwash or other antibacterial solutions. Ask your dentist to recommend an appropriate product for you.

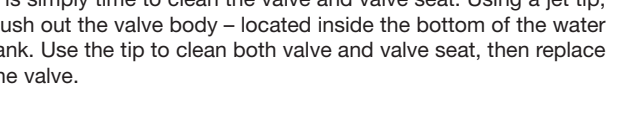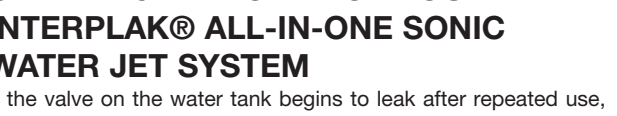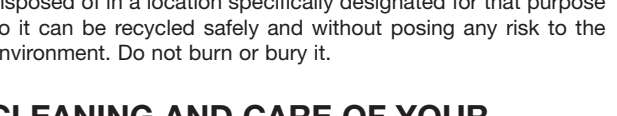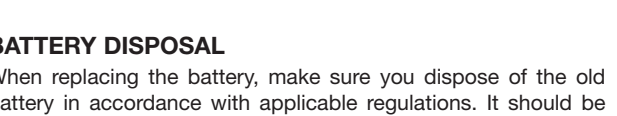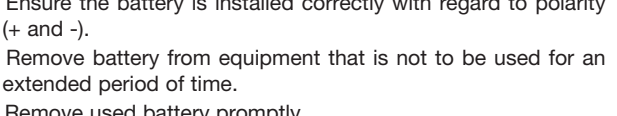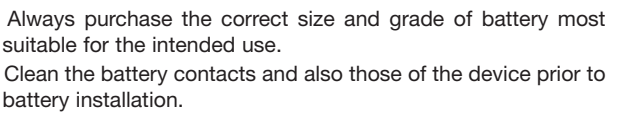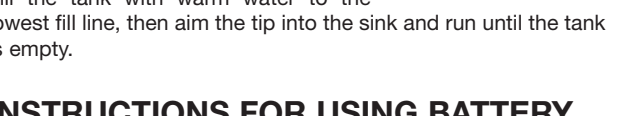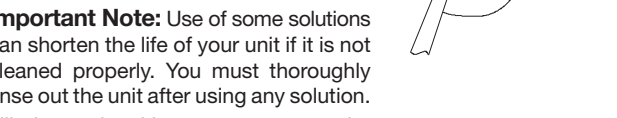TROUBLESHOOTING GUIDE

PROBLEM	CAUSE(S)	SOLUTION(S)
Sonic vibration not working	1. Dead battery 2. Dirty contact 3. Battery cap is loose 4. Sonic ON/OFF button on handle is in OFF position	1. Replace battery 2. Clean contact with cotton swab 3. Turn the battery cap to Closed position 4. Press Sonic ON/OFF button to ON
Tank leaks	The valve and valve seat in water tank need to be cleaned	Follow cleaning and care instructions on page 9
Inadequate pressure	1. Tank is not positioned properly 2. Tip is clogged with debris 3. Pressure control is on low setting	1. Reposition the tank and push down firmly 2. Try a different tip; clean or replace tip that is clogged 3. Adjust pressure control to higher setting
Main unit doesn't start	No power to the unit	Check the outlet and make sure it is functioning; check electrical reset buttons and wall switches
Still have questions? Call (02)8751-3633		

RECOMMENDED TIP REPLACEMENT INTERVALS

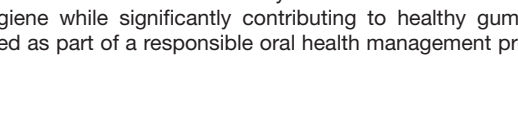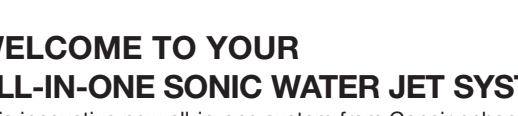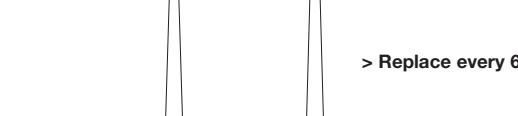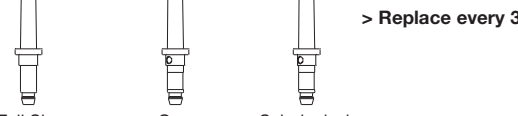THE INTERPLAK DIFFERENCE

The Interplak® All-in-One Sonic Water Jet System by Conair is designed to help you improve and maintain the highest level of oral hygiene quickly and easily, right at home. This integrated system gently and effectively cleans and flushes food debris from hard-to-reach places – around orthodontic appliances, fixed and removable bridgework, crowns, splints and implants – and provides valuable gum massage with every use.

BENEFITS OF SONIC TECHNOLOGY

Your Interplak® All-in-One Sonic Water Jet System by Conair is the only water jet with sonic vibrating action. All-in-one sonic technology provides enhanced brushing, flossing and massage, all in one system. The simultaneous sonic vibration and water pulsation greatly expand the performance of the dental water jet, offering a unique and complete cleaning experience and providing the ultimate in oral care.

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING

- Always unplug after use. This sonic water jet system should never be left unattended when plugged in.
- Do not place or store the unit where it can fall or be pulled into the tub, sink or toilet.
- Do not immerse unit in water or any other liquid.
- Do not reach for a unit that has fallen into water. Unplug immediately.**
- Do not reach into water under any circumstance. Have the unit checked by an Authorized Service Center after it has been immersed in water.
- Do not use while bathing.
- Do not handle with wet hands.
- The unit is designed for use with 110V.

9. Never operate the unit if it has a damaged power cord, if it is not working properly, if it has been dropped or damaged or if it has been immersed in liquids.

- Do not excessively bend or pull the power cord. Unplug the unit by grasping the plug – not the power cord.
- Keep the power cord away from heated surfaces.
- Close supervision is necessary when this product is used on or by children or disabled persons.
- Use this product only for its intended use as described in this manual or recommended by your dental professional.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product can produce pressures that can damage sensitive areas. See Recommended Technique for correct usage.

WELCOME TO YOUR ALL-IN-ONE SONIC WATER JET SYSTEM

This innovative new all-in-one system from Conair enhances oral hygiene while significantly contributing to healthy gums when used as part of a responsible oral health management program.

- Never block the air openings of the sonic water jet or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, and the like.
- Never drop or insert any object into any opening or into the hose.
- Fill reservoir with lukewarm water only unless otherwise specified by your dental professional. Do not overfill.
- Only use attachments recommended by Interplak.
- Never use while drowsy.
- Never operate where oxygen or aerosol sprays are being used.
- For household use only. Do not use outdoors.

SAVE THESE INSTRUCTIONS

MEDICAL WARNING For your safety:

- If you have had oral or periodontal surgery in the past 12 months, contact your dentist prior to use.
- When you start any new oral hygiene program, your gums may initially bleed slightly. Contact your dentist if excessive bleeding occurs or if the slight bleeding continues after two weeks of use.
- Do not use mouthwashes or rinses that contain iodine, bleach or tree oil.
- Do not use while wearing oral jewelry. Remove prior to use.
- Do not use if you have an open wound on your tongue or in your mouth.
- If your physician or cardiologist has advised you to receive antibiotic premedication before a dental procedure, you should consult your dentist before using this oral hygiene aid.
- Closely instruct and supervise children and individuals with special needs about proper use and prevent them from playing with this appliance.
- Do not remove the flexible tip on either the massage or the subgingival attachments. If it separates for any reason, discard both the tip and shaft and replace them with new attachments.
- Be sure that only the water jet tips touch your teeth, gums or tongue. Damage to those areas may occur if any other part of the handle touches them.
- Use only for oral hygiene cleaning. Do not use for any other purpose.
- As with any new oral hygiene program, individuals with heart problems should check with their physician prior to using the Interplak Sonic Water Jet System.

WARNING: Changes or modifications to this unit not expressly approved by Interplak will void the warranty. Servicing should always be referred to an authorized service center.

ONE YEAR LIMITED WARRANTY

This Conair Appliance is guaranteed for consumer use for one year. Please fill in and return the Warranty Registration as instructed or log into www.conair-asia.com to register. This guarantee is no longer valid in case of alteration or repair by any unauthorized persons and only covers defects under normal use from the date of purchase. If your product does not perform satisfactorily because of defects in materials or manufacture, it will be repaired or replaced under guarantee through the Conair Service Centre. Conair reserves the right to amend terms and specifications without prior notice. To obtain service under this warranty, please send in registration within one month of purchase through either one of the following ways.

ON-LINE REGISTRATION

Please log on to www.conair-asia.com.

OUT OF GUARANTEE

Call agents in your region as listed in this Instruction Booklet for details of repair service.

CUSTOMER SERVICE:

1F, No.417, Sec. 2, Tiding Blvd., NeiHu Dist., Taipei City 114, Taiwan (R.O.C.)
Conair Customer Hotline: (886) 2-8751-3633
E-mail: email_asia@conair.com
Website: www.conair-asia.com